

Model : AM4013MT Dino-Lite Premier

- 10X ~ 50X , 200X
- Enhanced 1.3 Megapixels
- With measurement and calibration
- Aluminum alloy case



Metal




Measurement and Calibration Capabilities



1.3 megapixel sensor



White Light

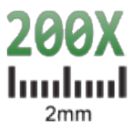

Product Description			
Product Number	AM4013MT	Magnification	10X ~ 50X , 200X
Warranty Period	1	Pixels	1.3 M pixels
Condition	New	Resolution	1280×1024
Packaging	Retail Box	Data Output	USB 2.0
Sensor	Enhanced Color CMOS	Video frame Rate	up to 30 frame-per-second
Illumination	8 white LED lights	Save File format	BMP, JPG, AVI
LED On/Off Function	Yes	Calibration Function	Yes
Body length	4 inches / 10.2 cm	Body Diameter	1.26 inches / 3.2cm max.
Cable length	6 feet / 180 cm	Unit	90(g)

Product Applications

The Dino-Lite Premier AM4013MT Premier digital handheld microscope features with enhanced 1.3MP sensor, providing a microscope that has high magnification at a nice distance with the ability to reduce glare and have advanced image quality.

The aerospace-grade machine tooled aluminum housing offers enhanced protection and durability even in the harshest working environments and the precision machined construction allows for smoother adjustment and operation and focusing. The stylish and durable polished aluminum body and performance upgrades make the AM4013MT the first choice for professionals.

The AM4013MT enhanced 1.3 megapixel sensor provides improved low light visibility, so you can view and see more detail when lighting is not perfect. The enhanced 1.3 megapixel sensor also provides very detailed images at its max resolution of 1280×1024 so you do not miss a thing.

The AM4013MT comes with the DinoCapture 2.0 software for PC's and DinoXcope software for MAC's. It allows you to take pictures, record videos, and also annotate on images. With professional models such as this, there is also measurement and measurement calibration tools to instantly gather or share important information of what you're looking at.

The Dino-Lite Premier AM4013MT digital handheld microscope is the perfect companion for people who require high magnification at a distance and the ability to reduce glare from highly reflective objects. Plug the USB to your computer and use the included software. The AM4013MT digital handheld microscope provides the ultimate digital microscope experience