

NEW PRODUCT

Model : AM4013ZT4 Dino-Lite Premier POLARIZER




- Handheld digital microscope with polarizer flyer
- 400x-470x magnification (non variable magnification model)
- Adjustable polarization
- USB connection
- Enhanced 1.3MP images
- Measurements and calibration



Skin without polarized



Skin with polarized

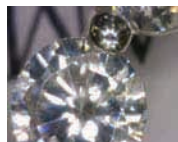

Jewelry without polarized



Jewelry with polarized



Turn Polarize off



Turn Polarize on

Product Description

Product Number	AM4013ZT4	Magnification	400x~470x
Warranty Period	1	Pixels	1.3M (SXGA)
Condition	New	Resolution	1280x1024
Packaging	Retail Box	Data Output	USB 2.0
Sensor	Enhanced Color CMOS	Video frame Rate	up to 30 frame-per-second
Illumination	built-in 8 white LED's	Save File format	BMP, JPG, AVI
LED On/Off Function	Yes	Calibration Function	Yes
Body length	4 inches / 10.2 cm	Body Diameter	1.26 inches / 3.2cm max.
Cable length	6 feet / 180 cm	Total package weight	280g +/- 20g

Product Applications

The Dino-Lite AM4013TZ4 digital handheld microscope has 400x ~ 470 x magnifications with adjustable polarization capabilities. Allows you to see fine details of surfaces by fully or partially removing the reflection of semi-transparent materials like; jewelry, skin, glass, or plastic. It is also deal for examining highly reflective materials, such as metals or solder joints. The AM4013ZT4 is part of the Dino-Lite Premier series that features enhanced 1.3 megapixel sensor for improved image quality.

Connect the Dino-Lite AM4013TZ4 to your computer via USB and use it with the included DinoCapture 2.0 software. Allowing you to take pictures, record video and also annotate right away. Having a Dino-Lite handheld digital microscope is a great tool to easily and instantly see magnified objects right on your computer screen.

Giving you ease of use and unique handheld microscopy experience.