




## Model : BL-ZW1 Dino-Lite Backlight Pad

- Best results with AM-413ZT, AM-412ZT or AD-413ZT
- Backlight stage with Polarizer
- Attachable target holder included
- Simple polarizer adjustment
- Built in illumination powered via USB or AC Adapter

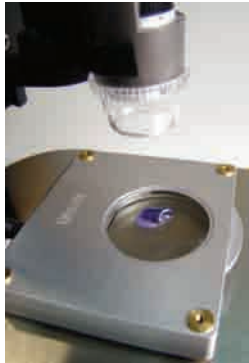

### Product Description

Product Number	BL-ZW1	Dimension	10cm x 9.62cmX2.05cm
Warranty Period	1	Type	Backlight Pad
Condition	New	Unit Weight	330 gm

### Product Applications

The BL-ZW1 is a back light illumination table with adjustable polarization for the Dino-Lite.

Please note that the polarizing feature has to be used with the Dino-Lite digital microscope with polarizing abilities; such as the AM413ZT or the AD413ZT handheld digital microscope. An optional accessory that can be attached to the BL-ZW1 is the BL-GC1 Gem Clip (The BL-GC1 Gem Clip is sold separately).

The polarizing ability of the BL-ZW1 can be used for many applications such as but not limited to geology (examining thin rock sections), metallurgy, and chemistry. For a business owner point of view, it can be used for jewelers for inspecting and presenting jewelry/ gems. The BL-ZW1 polarizing ability can also be useful for industrial applications; such as inspecting flexible circuit boards (FPC) for flaws and defects.

Power Source : USB or DC adaptor switchable

Operating Voltage : DC 5V± 0.5V

Light Intensity : over 5000 cd /m<sup>2</sup> Freely rotatable polarizer works with handheld digital microscope AM413ZT or handheld digital microscope AD413ZT

Housing Material : Aluminum Alloy